


EMTR Series TRINOCULAR STEREO BODY


MEIJI TECHNO EMTR Series is not only a two magnification type stereo microscope body with turret mounted objective pairs, but also provides a photo tube for use in television and photographic applications.

When the photo prism is drawn into the optical path, 100% of the light is directed to the photo tube and the image can still be examined visually through the right-hand eyepiece.

- Crisp, erect image with flat, widefield of view and long working distance.
- Paired widefield objectives in a turret mount for rapid selection of magnification.
- Binocular tubes with an inclination of 45° for ease of observation.
- Trinocular body can be rotated a full 360° and locked in any position.
- Interpupillary distance adjustment from 54mm to 75mm.
- Dioptric adjustment is possible on the left eyepiece tube.
- Vertical photo tube is provided for easy compatibility with CCTV or photographic system.

IT PROVIDES COMPLETE VERSATILITY AND MOUNTING COMPATIBILITY WITH THE PRESENT LINEUP OF MEIJITECHNO EMT SERIES STEREO MICROSCOPE ACCESSORIES.

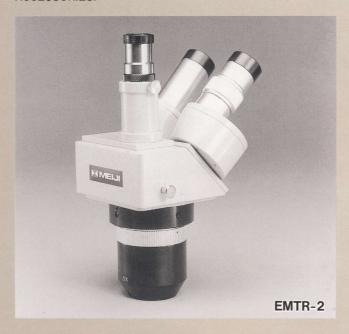

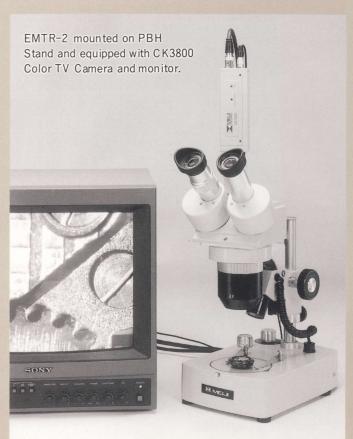

ORDERING INFORMATION

MODEL	OBJECTIVES	EYEPIECES
EMTR-1	IX and 2X	SWF10X, paired
EMTR-2	IX and 3X	SWF10X, paired
EMTR-3	2X and 4X	SWF10X, paired
EMTR-4	IX and 4X	SWF10X, paired

EVEDIECES & EVECUIEI DO

ETEPIEC	ES & ETESHIELDS
MA501	SWF5X Super widefield, paired, Field No. 26
MA502	SWF10X Super widefield, paired, Field No. 23
MA519	SWF10X Focusing eyepiece
MA520	SWF12.5X Super widefield, paired, Field No. 20
MA503	SWF15X Super widefield, paired, Field No. 15.4
MA535	HSWF15X Super widefield, high eyepoint, paired,
	Field No. 16
MA504	SWF20X Super widefield, paired, Field No. 11.5

SWF30X Super widefield, paired, Field No. 7.7 MA521

MA600 New Eyeshield, paired

SUPPLEMENTARY LENSES

MA518	Supplementary lens 0.5X for EMTR-2, EMTR-3 & EMTR-4
	(MA515 Adapter is needed to use MA518 on EMTR-2)
MA527	Supplementary lens 0.75X for EMTR-2, EMTR-3 & EMTR-4
	(MA515 Adapter is needed to use MA527 on EMTR-2)
MA513	Supplementary lens 1.5X for EMTR-2, EMTR-3 & EMTR-4
	(MA515 Adapter is needed to use MA513 on EMTR-2)
MA514	Supplementary lens 2X for EMTR-2, EMTR-3 & EMTR-4
	(MA515 Adapter is needed to use MA514 on FMTR-2)

MA515 Adapter for EMTR-2 (use to attach MA518, MA527, MA513 & MA514 on EMTR-2)

EMTR Series OPTICAL DATA

	OBJEC- TIVE	EYEPIECES	SWF10X		SWF15X		SWF20X		SWF30X		Working
MODEL		SUPPLEMEN- TARY LENS	Total mag.	Field Diam (mm)	Total mag.	Field Diam (mm)	Total mag.	Field Diam (mm)	Total mag.	Field Diam (mm)	
EMTR-1	1X, 2X	-	10X, 20X	23.0, 11.5	15X, 30X	15.3, 7.7	20X, 40X	11.5, 5.8	30X, 60X	7.6, 3.8	108
			10X, 30X	23.0, 7.7	15X, 45X	15.3, 5.1	20X, 60X	11.5, 3.8	30X, 90X	7.6, 2.5	81
EMTR-2	1X, 3X	0.5X	5X, I5X	46.0, 15.3	7.5X, 22.5X	30.6, 10.2	10X, 30X	23.0, 7.6	15X, 45X	15.3, 5.1	96
		2.0X	20X, 60X	11.5, 3.8	30X, 90X	7.6, 2.5	40X, I20X	5.7, 1.9	60X, 180X	3.8, 1.2	21
	2X, 4X	_	20X, 40X	11.5, 5.7	30X, 60X	7.6, 3.8	40X, 80X	5.7, 2.8	60X, I20X	3.8, 1.9	63
EMTR-3		0.5X	10X, 20X	23.0, 11.5	15X, 30X	15.3, 7.6	20X, 40X	11.5, 5.7	30X, 60X	7.6, 3.8	96
		2.0X	40X, 80X	5.7, 2.8	60X, I20X	3.8, 1.9	80X, 160X	2.8, 1.4	120X,240X	1.9, 0.9	21
	IX, 4X	_	10X, 40X	23.0,5.7	15X, 60X	15.3,3.8	20X, 80X	11.5,2.8	30X, 120X	7.6,1.9	63
EMTR-4		0.5X	5X, 20X	46.0,11,5	7.5X, 30X	30.8,7.6	10X, 40X	23.0,5.7	15X, 60X	15.4.3.8	96
		2.0X	20X, 80X	11.5,2.8	30X, I20X	7.6,1.9	40X, 160X	5.7,1.4	60X, 240X	3.8,0.9	21

TV EQUIPMENT

CCD TV CAMERAS

A separate "C" mount adapter is required to connect the CCD TV camera to the microscope.

MK300N	CCD B/W TV camera, I/2" CCD, 570 TV lines of horizontal resolution,
	NTSC system, with 115V-DC12 adapter

CCD B/W TV camera, I/2" CCD, 570 TV lines of horizontal resolution, MK300P PAL system, with 220V-DC12 adapter

CCD Color TV camera, I/2" CCD, 450 TV lines of horizontal resolution, CK3800N NTSC system, with 115V-DC12V adapter

CCD Color TV camera, 1/2" CCD, 450 TV lines of horizontal resolution, PAL system, with 220V-DCI2V adapter

A detailed catalog of the CK3800 CCD Color TV camera is available upon request.

"C" MOUNT ADAPTERS

MA151/5N	"C" mount adapter to connect EMTR-1, EMTR-2, EMTR-3, EMTR-4 to TV camera.
----------	---

"C" mount adapter with 1.0X lens (for use with existing photo tubes of MA151/35/15 EMTR-1, EMTR-2, EMTR-3, EMTR-4)

"C" mount adapter with 0.7X lens (for use with existing photo tubes of MA151/35/20 EMTR-1, EMTR-2, EMTR-3, EMTR-4)

"C" mount adapter with 2.5X lens (for use with existing photo tubes of MA151/35/25 EMTR-1, EMTR-2, EMTR-3, EMTR-4)

Agent: MEIJI TECHNO AMERICA



5895 Rue Ferrari, San Jose, CA 95138, USA Tel: 408-226-3454 Fax: 408-226-0900 Toll Free: 1-800-832-0060

E-mail: info@meijitechno.com http://www.meijitechno.com



CK3800 1/2"CCD COLOR TV CAMERA

Comprises:	
CK3800 color TV camera	
AC Adapter ·····	I рс
BNC Cable, 3M long ·····	l pc
S Cable, 3M long ·····	I рс

MEIJI TECHNO CO.,LTD.

322-1, Chikumazawa Miyoshi-machi, Iruma-gun Saitama 354-0043, Japan

E-mail: meiji@meijitechno.co.jp Phone: (0)49-259-0111 Fax: (0)49-259-0113