MEJ L3000 SERIES METALLURGICAL MICROSCOPES

Infinity Optical System · Reflected & Transmitted Brightfield Illumination

FEATURES VIEWING HEAD:

The ML8000 Series Metallurgical Microscope is equipped with either an erect image binocular or regular trinocular viewing head inclined at 30° and rotatable through 360° for comfortable viewing. Each eyetube is provided with a graduated dioptre adjustment. Interpupillary distance is adjustable between 52mm and 72mm.

EYEPIECES:

Paired DIN KHW10X super wide 20mm field of view evepieces.

OBJECTIVE CHANGER:

Smooth operating quintuple ball-bearing nosepiece.

OBJECTIVES:

Infinity corrected, S. Planachromat M4X/0.10, M10X/ 0.25, M20X/0.40, M40X/0.65


Vertical full Koehler illuminator with 6V 30W halogen lamp for brightfield and polarization technique, with infinity lens. Electronic transformer with intensity control supplied in a separate case. Blue clear, green clear, ND8 Neutral density and polarizing fil-

Transmitted Koehler illuminator with 6V 30W halogen lamp with electronic transformer and intensity control built in the microscope base.

Large scan 4" x 4" ball-bearing mechanical stage with 94.5mm diameter clear glass stage plate, with drop-down coaxial axis controls with X/Y movement.

Rack and pinion focusable substage condenser with iris diaphragm and filter tray, including clear blue

Coaxial low position coarse and fine focus controls, graduated to 2 micron per division, are fitted with tension adjustment and safety autostop. The focusing range is 35mm. The movement runs in extra wide

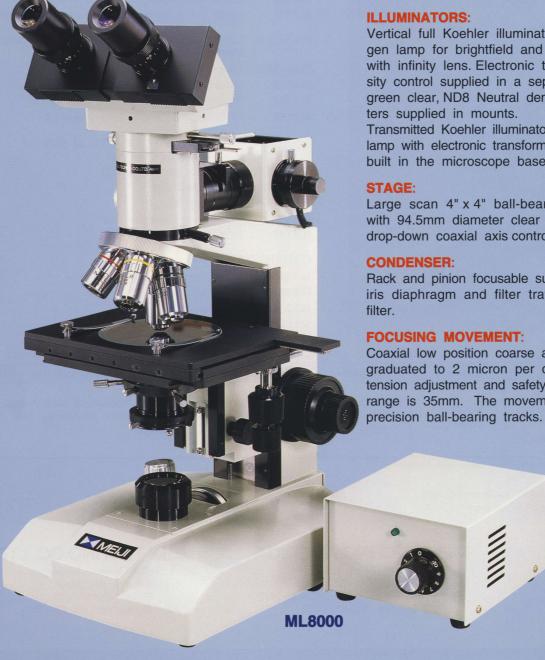

VIEWING HEADS

MA401	Binocular head, erect image (included	d in
	Model ML8000)	

MA603 Trinocular head (included in Model ML8100)

MA495 Trinocular head, erect image

MA403 Monocular head

EYEPIECES (DIN)

MA407	KHW10X Eyepiece, compensating, FN 20mm	
	(included)	

MA406 HWF10X Eyepiece, High eyepoint, FN 18mm
 MA408 HWF15X Eyepiece, High eyepoint, FN 12.2mm
 MA409 HWF20X Eyepiece, High eyepoint, FN 9mm



OBJECTIVES (Infinity corrected)

MA336	M4X/0.10 S. Plan, W.D. 20mm (included)
MA330	M5X/0.10 S. Plan, W.D. 20mm
MA337	M10X/0.25 S. Plan, W.D. 9.4mm (Included)
MA338	M20X/0.40 S. Plan, W.D. 5.28mm
MA339	M40X/0.65 S. Plan, W.D. 0.81mm (Included)
MA344	M50X/0.70 S. Plan, W.D. 0.36mm
MA349	LWD M50X/0.50 S. Plan, W.D. 6.0mm
MA329	M60X/0.85 S. Plan, W.D. 0.38mm
MA328	M100X/0.75 S. Plan, dry, W.D. 0.27mm

ACCESSORIES FOR MEASUREMENT

MA415	Filar micrometer	eyepiece	10X	with	inclined
	monocular head				

MA255 Eyepiece micrometer, 10mm into 100 parts

for KHW10X

MA256 Eyepiece micrometer, 5mm into 100 parts

for KHW10X

MA284 Eyepiece micrometer, crossline reticle for

KHW10X

MA285 Stage micrometer, 1mm into 100 units

MA292 Metal stage micrometer, 1mm into 100 units

FILTERS 15mm diam. in mounts for use with

Vertical Illuminator

MA475/05 Blue clear filter (included)
MA476 Green clear filter (included)
MA497 Polarizing filter (included)

MA198 ND8 Neutral density filter (included)

FILTERS 29.8mm diam. for use in tray of

Substage Condenser

MA289/05 Blue clear filter (included)

MA290/05 Green clear filter
MA291 Polarizing filter

MA196 ND2 Neutral density filter MA197 ND8 Neutral density filter

OTHER ITEMS

MA326 Spare lamp 6V 30W halogen

MA327 Fuse 3A MA701 Dust cover

MA133S Wooden cabinet with lock and key

ORDERING INFORMATION FOR ML8000 SERIES

XYour orders should always indicate your voltage for use.

Model	Body	Eyepiece	Objectives	Stage	Illuminator	
ML8000	Binocular (Erect image)	KHW10X 20mm f.o.v	S. Plan M4X, M10X M20X, M40X	Mechanical Stage, 4"x4", with clear glass Plate	Vertical Koehler, 6V 30W Halogen & Transmitted Koehler, 6V 30W	
ML8100	Trinocular				Halogen	

Supplied in styrofoam case with dust cover and manual.

SIZE AND WEIGHT:

 $\begin{array}{l} ML8000:227mm~(W)~x~415mm~(D)~x~472mm~(H),~6.8Kg\\ ML8100:227mm~(W)~x~415mm~(D)~x~500mm~(H),~6.9Kg \end{array}$

Agent:



MEIJI TECHNO AMERICA

5895 Rue Ferrari, San Jose, CA 95138, USA Tel: 408-226-3454 Fax: 408-226-0900

Toll Free: 1-800-832-0060 E-mail: info@meijitechno.com http://www.meijitechno.com



MEIJI TECHNO CO.,LTD.

322-1,Chikumazawa Miyoshi-machi, Iruma-gun Saitama 354-0043, Japan E-mail: meiji@meijitechno.co.jp

Phone: (0)49-259-0111 Fax: (0)49-259-0113