




HIGHLIGHTS

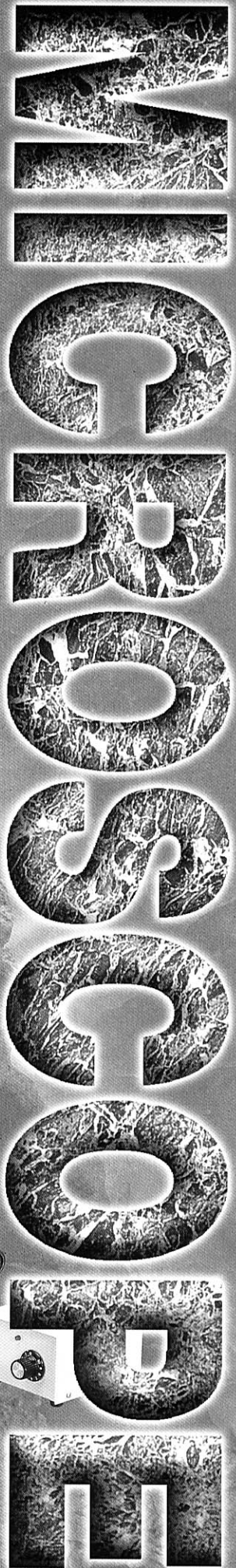

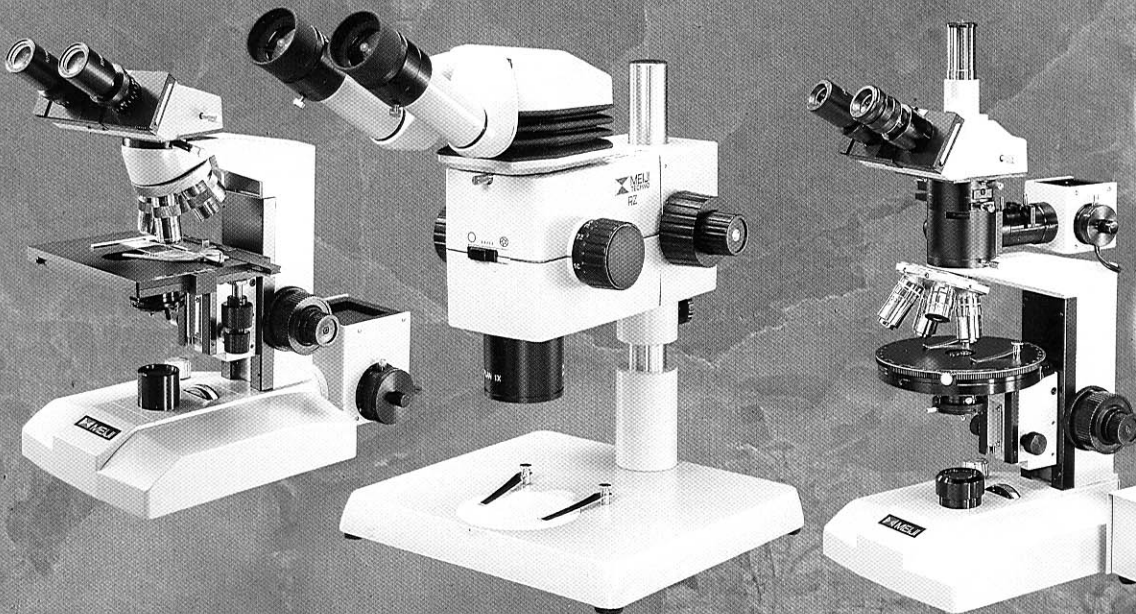


ADVANCED STEREOS

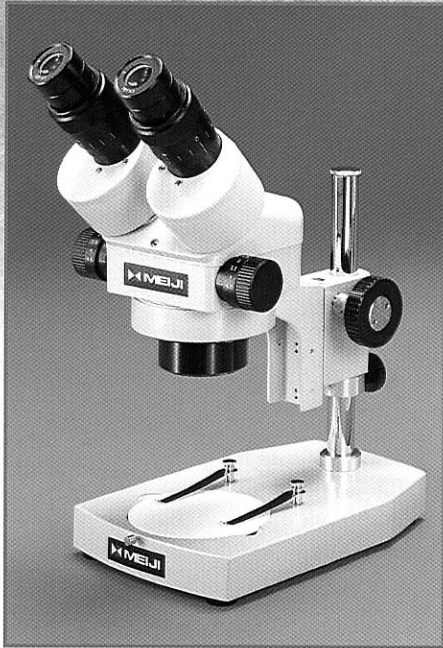


EMZ-5 Zoom Stereo Microscope

Stereo zoom body offering 0.7X - 4.5X zoom range, and having unequalled image quality and resolution, and with the advantage of side-mounted zoom controls.

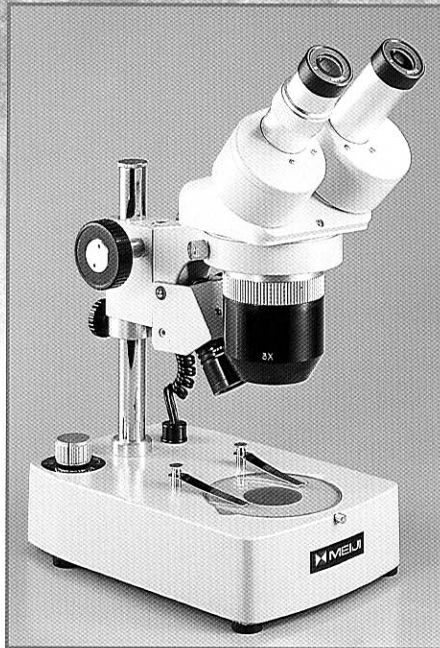


EMT Stereo Microscope

Dual magnification stereo microscopes with turret-mounted objective pairs. Total magnification ranges from 10X - 20X to 120X - 240X with varied and optional lenses.

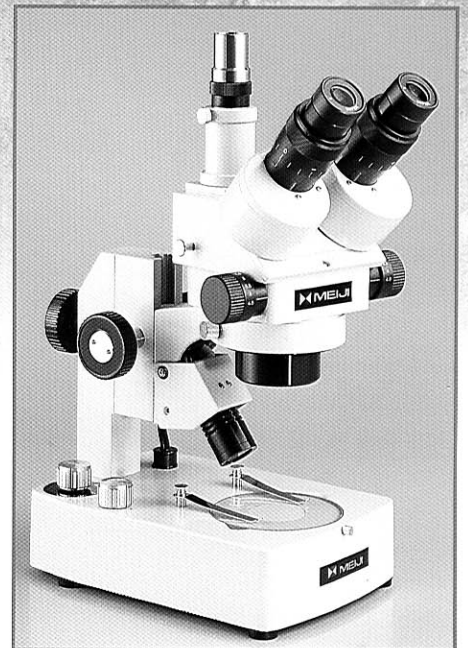


EMZ-5TR Trinocular Zoom Stereo Microscope

Zoom stereo trinocular body, optically the same as EMZ-5, but also offering full photographic and video capability.

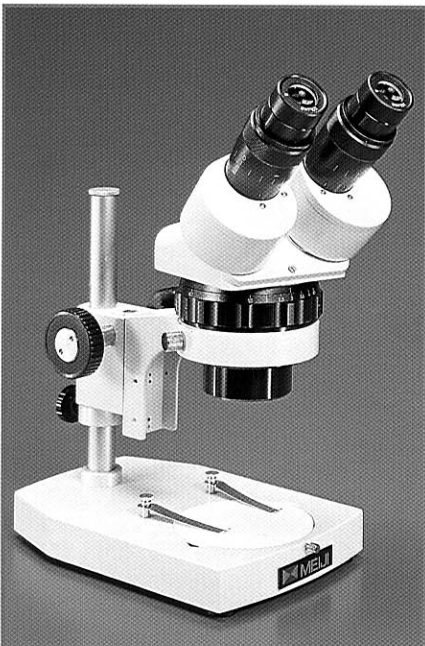


EMZ-2 Zoom Stereo Microscope

Stereo zoom body offering 0.7X - 4.5X zoom range, and having unequalled image quality and resolution with center zoom controls.

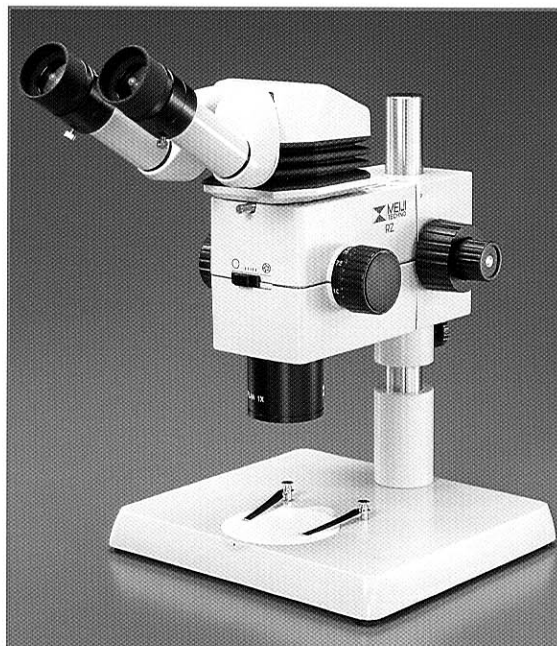

RZ Zoom Stereo Microscope

C.M.O. Type Zoom Stereo. Available in either ergonomic type or 45° inclined viewing head. Magnification from 3.75X to 300X through optional objectives and eyepieces.

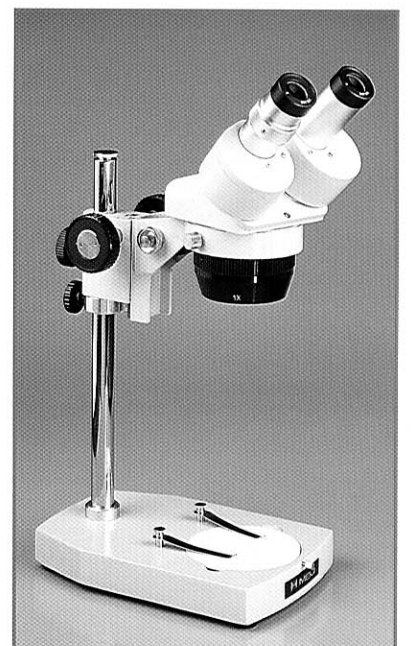


EMX Stereo Microscope, Long working distance

Long working distance stereo microscope with turret-mounted objective pair. Magnification ranges from 5X - 10X to 10X - 20X with optional eyepieces.




EDUCATIONAL MICROSCOPES

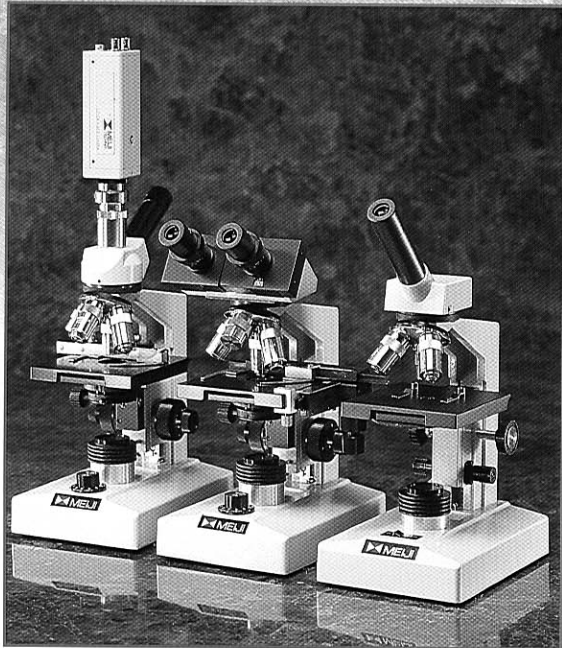


TM-400 Teaching Microscope

Advanced teaching microscopes with DIN optics, offering in monocular and binocular versions. Now also available with coaxial coarse and fine focusing.

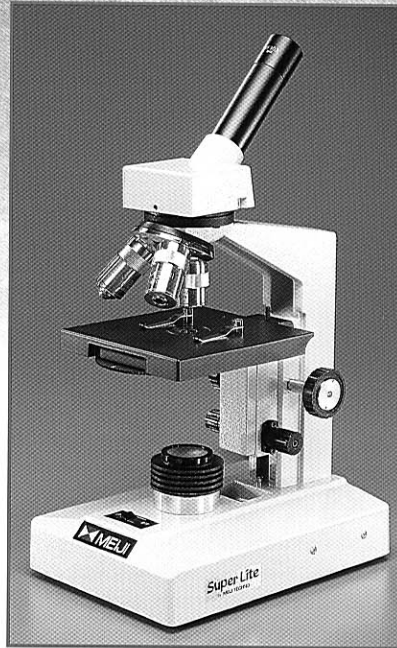


TM400F Teaching Microscope with cool "Super Lite"

Advanced teaching microscope with in-base "Super Lite" cool and long life fluorescent illuminator for 115V with 3-wire cord for use in USA.

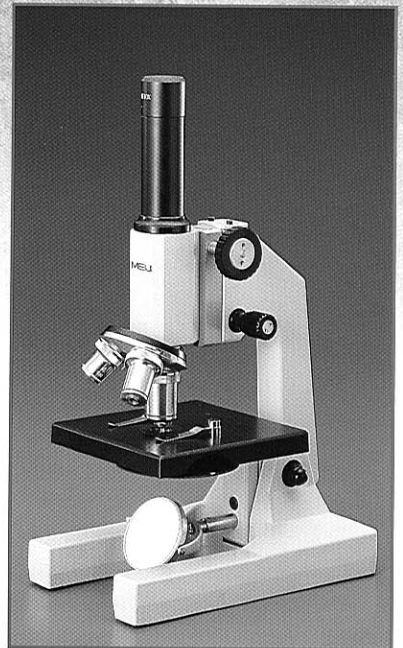

TM-200 Teaching Microscope

Most popular teaching microscope. Rugged, "Student-proof" design at an affordable price.


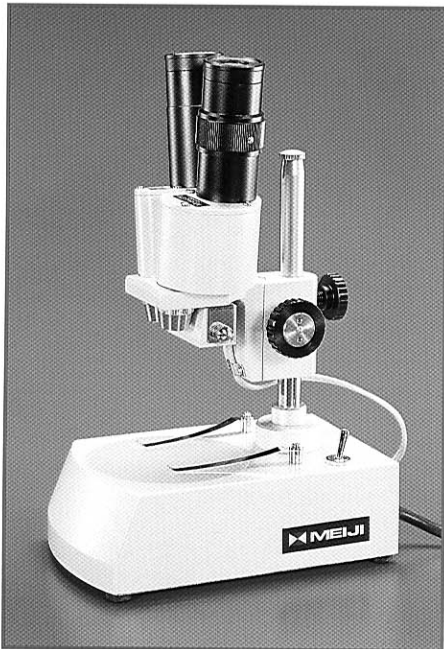


ST Stereo Microscope

An inexpensive stereo microscope with either fixed or interchangeable objective pairs. Magnifications available from 10X to 40X with optional eyepieces.

SK Stereo Microscope

Widefield non-modular stereo microscope with inclined head, and with objective pair in turret. Models with incident and transmitted lights are available.

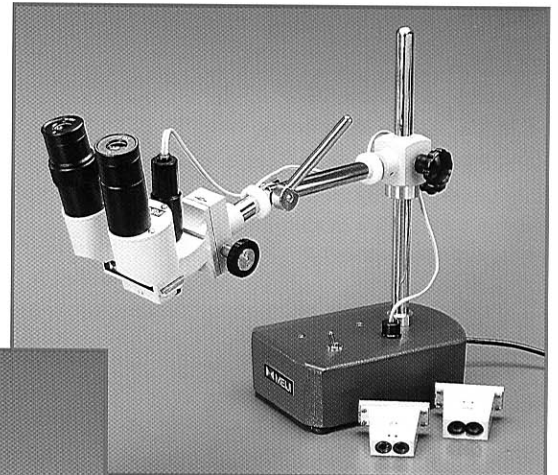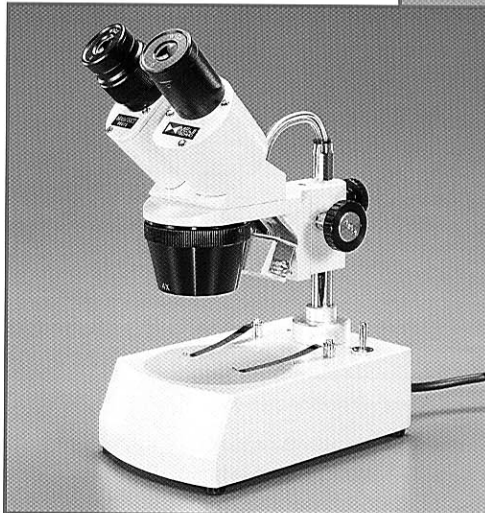


BM Stereo Microscope, long-arm stand

Stereo microscope mounted on heavy, long-arm stand with intense incident illuminator - 6V 2A. Power supply is built into the base.




COMPOUND HIGH-POWER MICROSCOPES

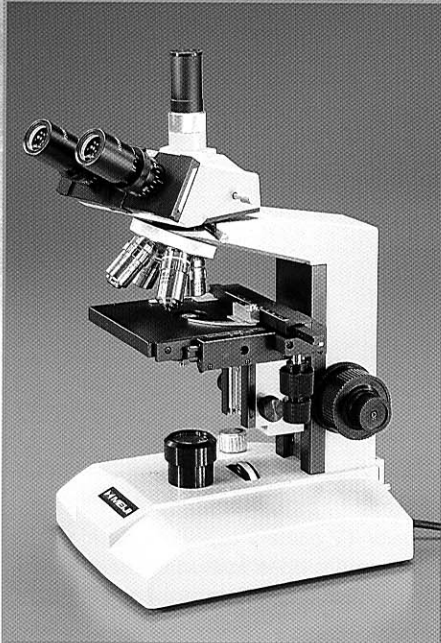


ML2000 Biological Microscope

Biological microscopes with Koehler type illuminator available with brightfield and Phase contrast optics. Either binocular or trinocular for full photographic and video applications.

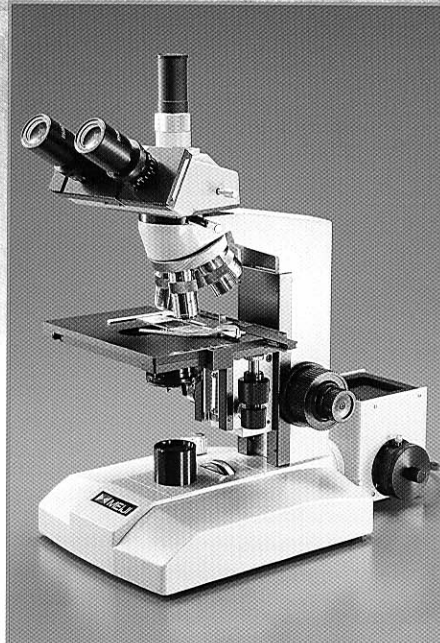


ML5000 Biological Microscope

A full size biological microscope with Koehler type Halogen Illuminator and a wide range of optional objective sets. Either binocular or trinocular for full photographic and video applications.

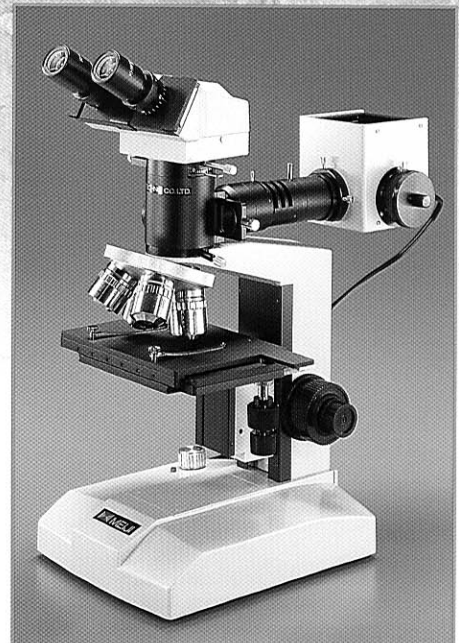


ML7500 Metallurgical Microscope

A Brightfield-Darkfield microscope with 12V 50W Halogen vertical illuminator. The newly designed optics offer superior image quality in both brightfield and darkfield settings. Full photo-graphic and video accessories are available.

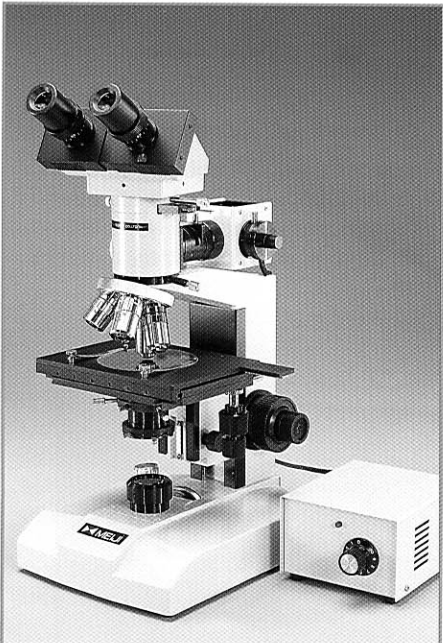


ML8000 Metallurgical Microscope

A brightfield metallurgical microscope with Koehler type vertical illuminator. Offered with a full set of filters including simple polar-izers. Either binocular or trinocular for photographic and video applications.

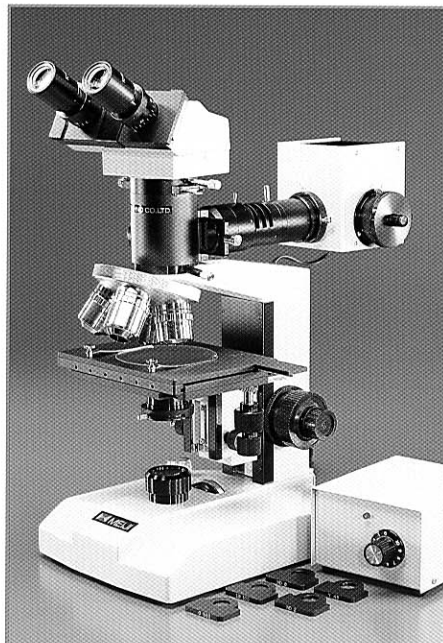

ML8500 Universal Microscope

A Brightfield and Darkfield microscope with both 12V 50W Halogen vertical illuminator and 6V 20W Halogen transmitted illuminator, and with substage condenser. Full photo-graphic and video accessories are available.

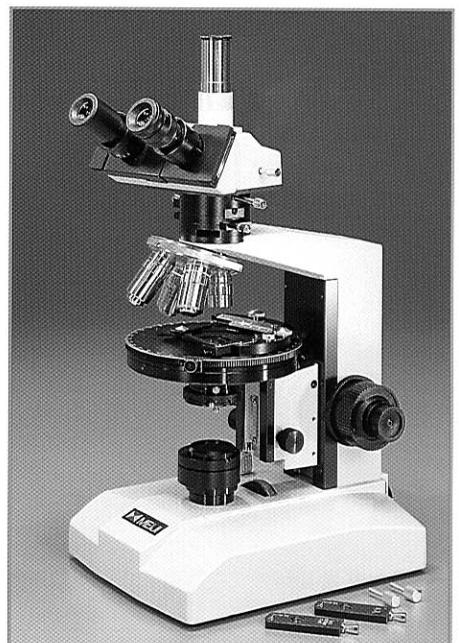


ML9000 Polarizing Microscopes

A full size polarizing microscope with centerable objectives and ball-bearing rotating stage. The Bertrand lens gives extra large and bright interference figures in both binocular and trinocular models.




SYSTEM-ACCESSORIES

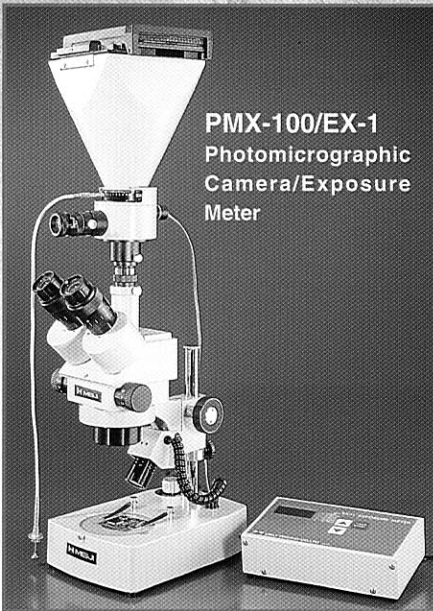


PMX-100/EX-1
Photomicrographic
Camera/Exposure
Meter

A large format camera body with manual shutter and viewing eyepiece, which accepts a variety of film holders. When using EX-1 semi-automatic exposure meter with PMX-100 accurate exposure time can be determined and set.

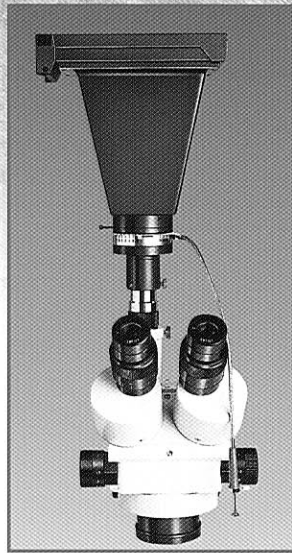

PMS-810 Polaroid
Photomicrographic Camera

Polaroid 3 1/4 x 4 1/4 camera with speed shutter and built-in 8X photo eyepiece.


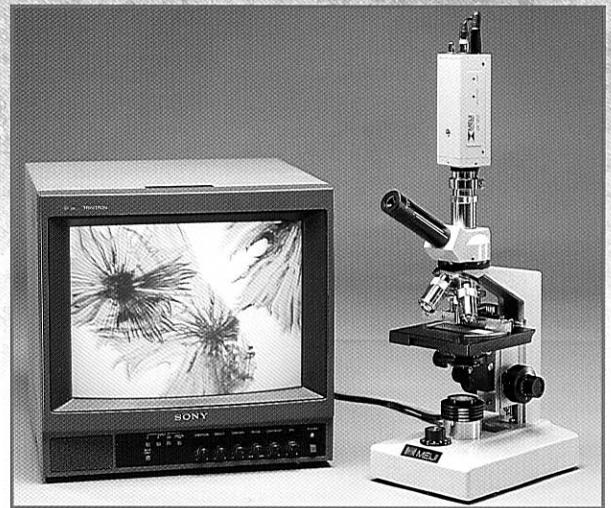


CK-3800 CCD Color TV Camera

A compact CCD Color TV camera especially designed for use with microscopes. Includes built-in automatic white balance control, and manual color adjusting controls.


FL180
Fiber Optic Illuminator

Fiber optic illuminator available with single, dual, or annular ring light guides. 150W of cold controllable light.

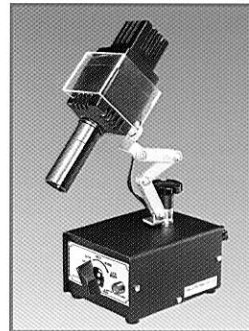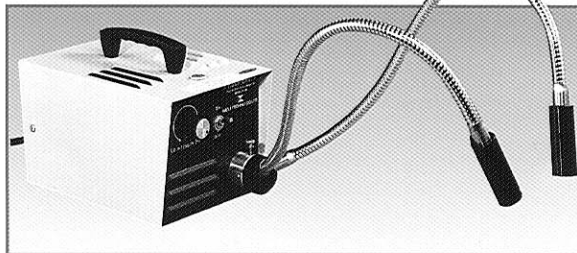


MA264
Versatile Illuminator

Intense illuminator 6.5V, 2.75A with prefocused 3-lens condensing system, mounted on detachable universal arm.

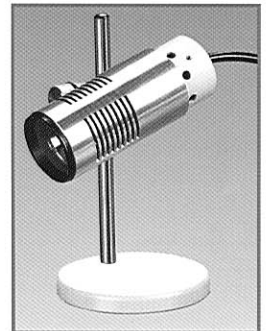


MA263 Free Stand
Illuminator

A high intensity illuminator with miniature flood bulb 30W, mounted on adjustable stand.

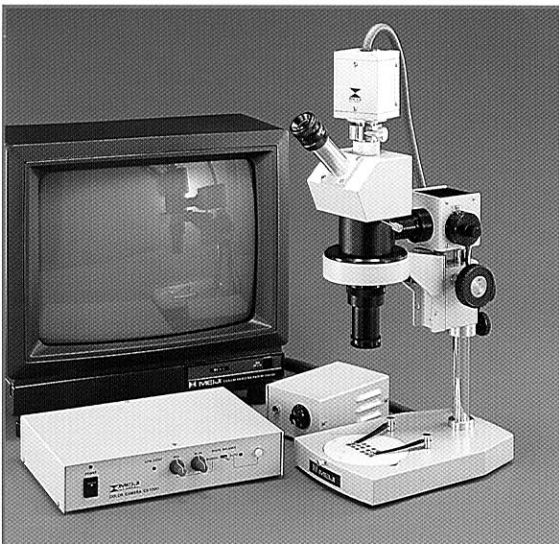

UNIMAC Zoom Macroscope

A fully modular universal zoom macroscope that lends itself to all applications including photomicrography and video display.

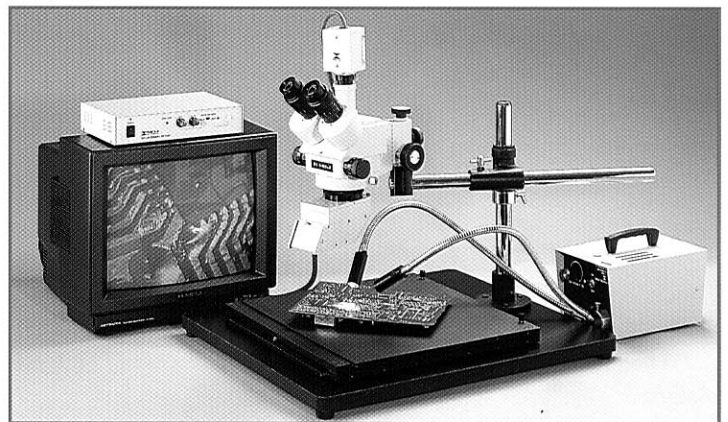

SMD Inspection System

Inspection system for inspection of mounted Flat package, LSI, chips and printed circuit board. Quick changeover between vertical and oblique observation can be made for viewing objects from all angles.




SYSTEM-ACCESSORIES




MA150/60 Camera attachment with finder

Camera adapter with intermediate viewing eyepiece for 35mm S.L.R. cameras.

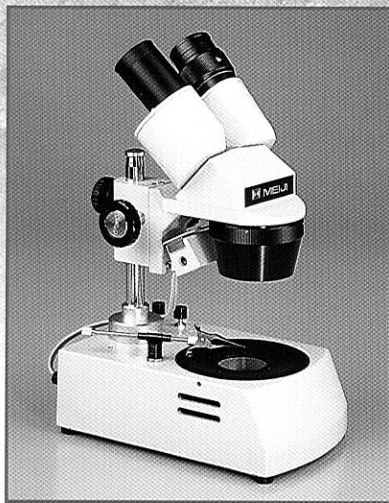


GS-2000 GEM Stereo Microscope

Stereo microscope, designed for examination of gem stones, with built-in darkfield illuminator and incident illuminator.

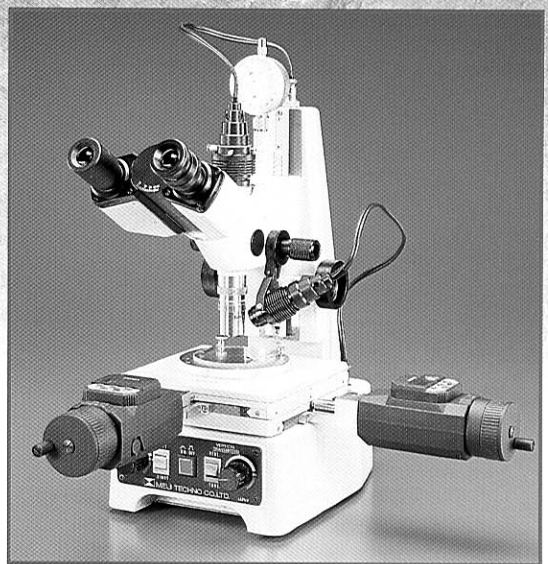

MT Measuring Microscope

Measuring microscope using infinity corrected optics and 30 degree inclined viewing head. A wide variety of accessories, including photomicrographic and video display, are available.


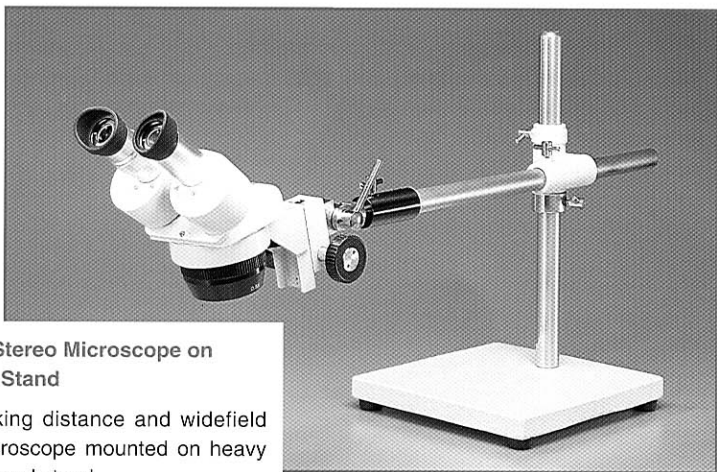


EMX-U Stereo Microscope on Universal Stand

Long working distance and widefield stereo microscope mounted on heavy base universal stand.

MA273DIN
Zernike Phase Contrast Set

A Phase contrast condenser with adjustable Phase annuli. Available with standard achromat or optional semi-Plan Phase objectives.

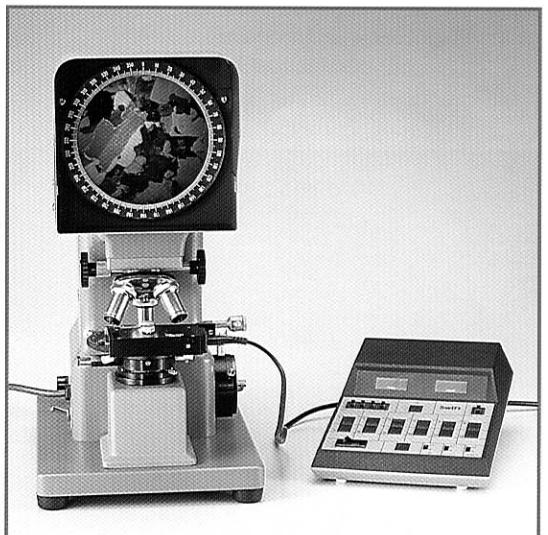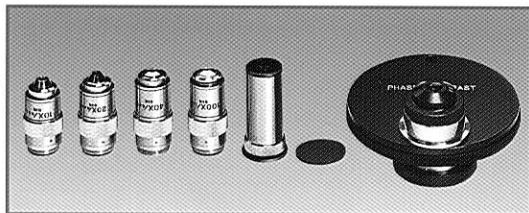

PV Polarizing Micro-projector

Polarizing microprojector optionally usable with automatic point counting equipment.

Agent



MEIJI TECHNO CO.,LTD.

6, Oi-670, Oi-machi, Iruma-gun
Saitama 356-0053, Japan
Phone : 492-67-0911

Fax : 492-69-0691, 492-69-0692

MEIJI TECHNO AMERICA

2186 Bering Drive
San Jose, CA., 95131, USA
Phone : 408-428-9654 Fax : 408-428-0472
Toll free : 800-832-0060